

Portland State University

PDXScholar

Student Research Symposium

Student Research Symposium 2021

Quantifying the motility of *Vibrio alginolyticus* after simulated microgravity using digital holographic analysis

Jacqueline Acres
Portland State University

Follow this and additional works at: <https://pdxscholar.library.pdx.edu/studentsymposium>

Let us know how access to this document benefits you.

Acres, Jacqueline, "Quantifying the motility of *Vibrio alginolyticus* after simulated microgravity using digital holographic analysis" (2021). *Student Research Symposium*. 5.
<https://pdxscholar.library.pdx.edu/studentsymposium/2021/Posters/5>

This Poster is brought to you for free and open access. It has been accepted for inclusion in Student Research Symposium by an authorized administrator of PDXScholar. Please contact us if we can make this document more accessible: pdxscholar@pdx.edu.


Quantifying motility of *Vibrio alginolyticus* after simulated microgravity using digital holographic microscopy

Jacqueline M Acres
Jay Nadeau
Portland State University

Introduction

- Microgravity is the environment on the space station.
- It's very difficult to conduct experiments on the space station, so we simulate microgravity on Earth using special devices. The ones used in this study were high aspect ratio vessels, a type of rotating wall vessel.
- Simulated microgravity mimics specific conditions:
 - Low fluid shear
 - Lack of sedimentation
- Motility gene expression after microgravity and simulated microgravity isn't consistent between spaceflight and simulated microgravity devices. To better understand this, a good first step is to quantify motility phenotypes.
- Digital holographic microscopy (DHM) is a non-invasive technique that allows us to image microbes in a 3D environment over time.

Methods

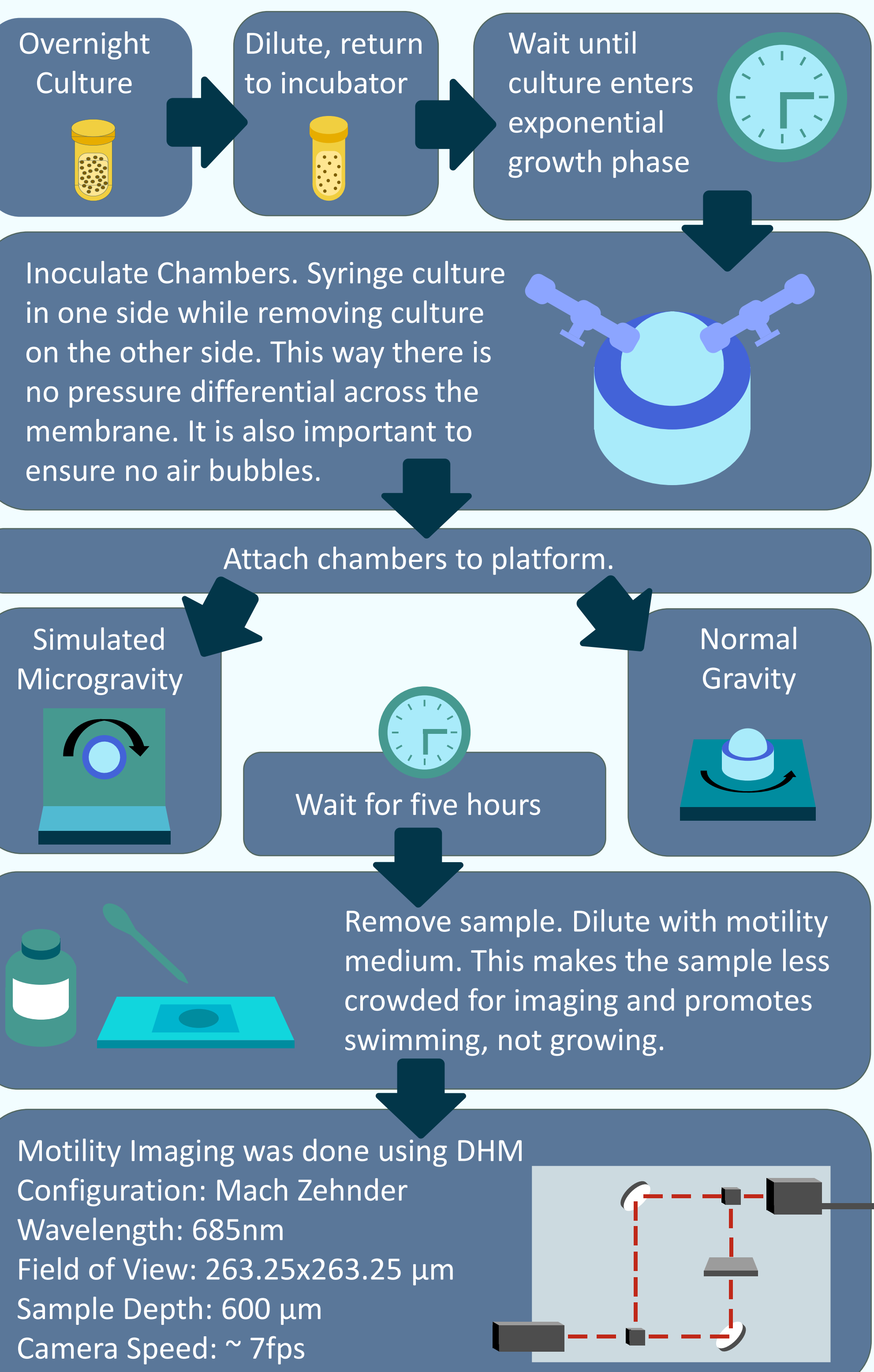


Reconstruct holograms. Concatenate z-planes into a hyperstack to scroll through time and z. Tracking was done manually. There is currently no high-throughput method to process DHM data.

How fast these bacteria move does NOT change after putting them in a weightless environment on Earth


Results / Analysis

1. Comparison between 2D and 3D traces of a microbe after normal gravity:

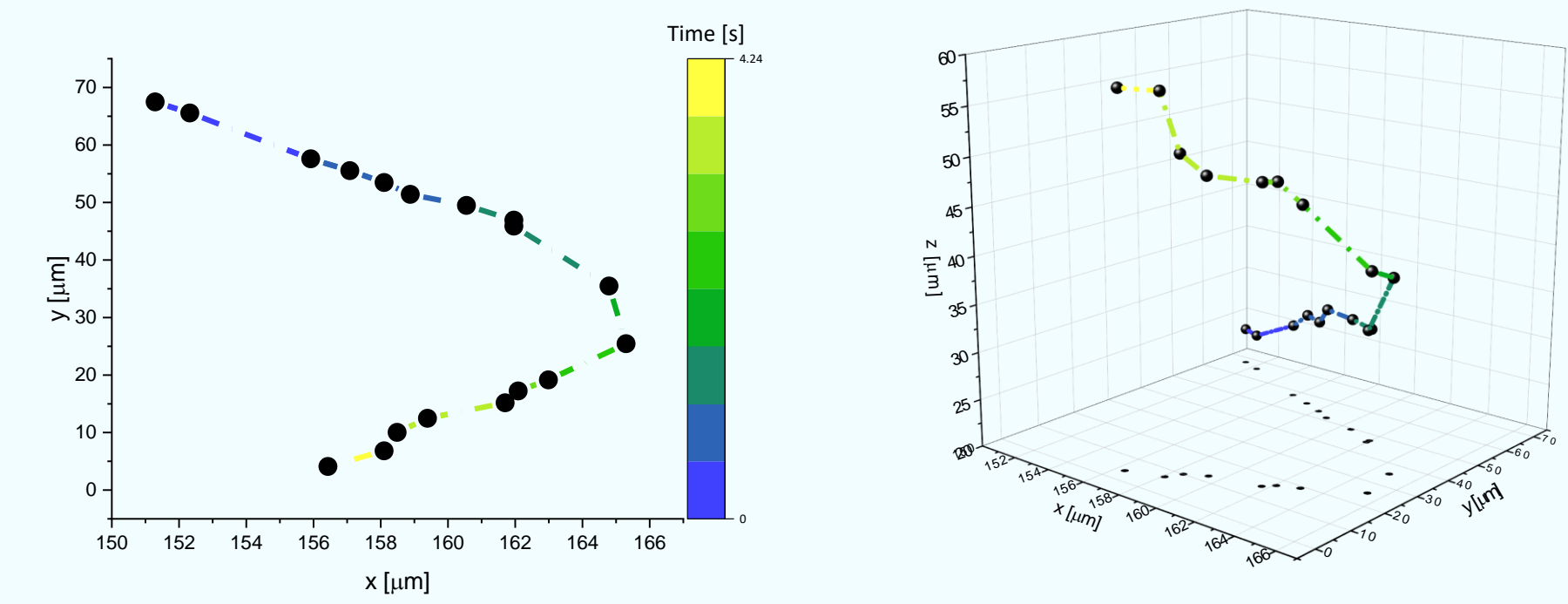


• 2D Traces can be misleading over organism movement. A 3D trace shows a clearer picture.

2. Run-reverse-flick motility of *Vibrio alginolyticus* observed via DHM

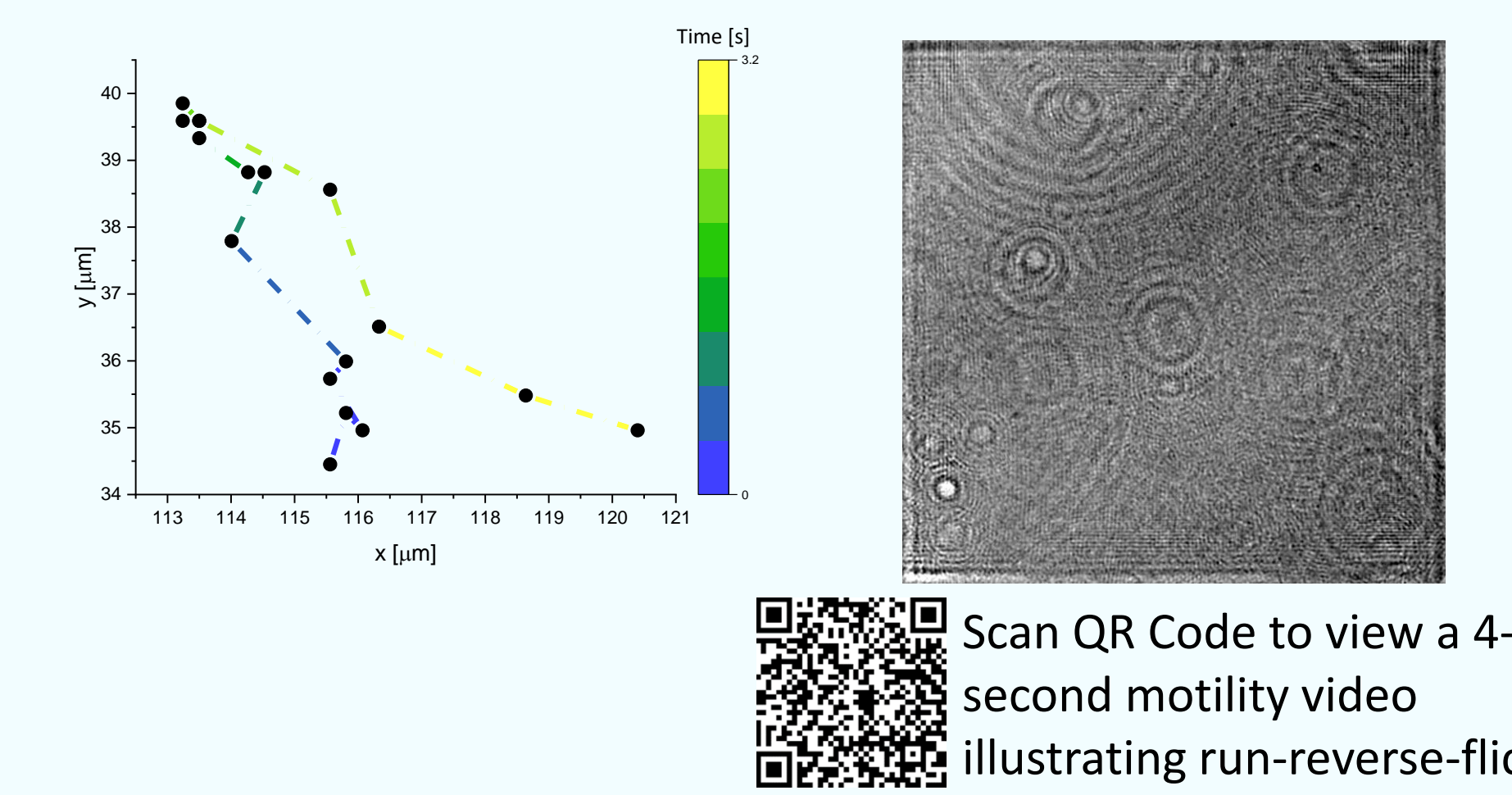


3. Impact of swimming speed a surface

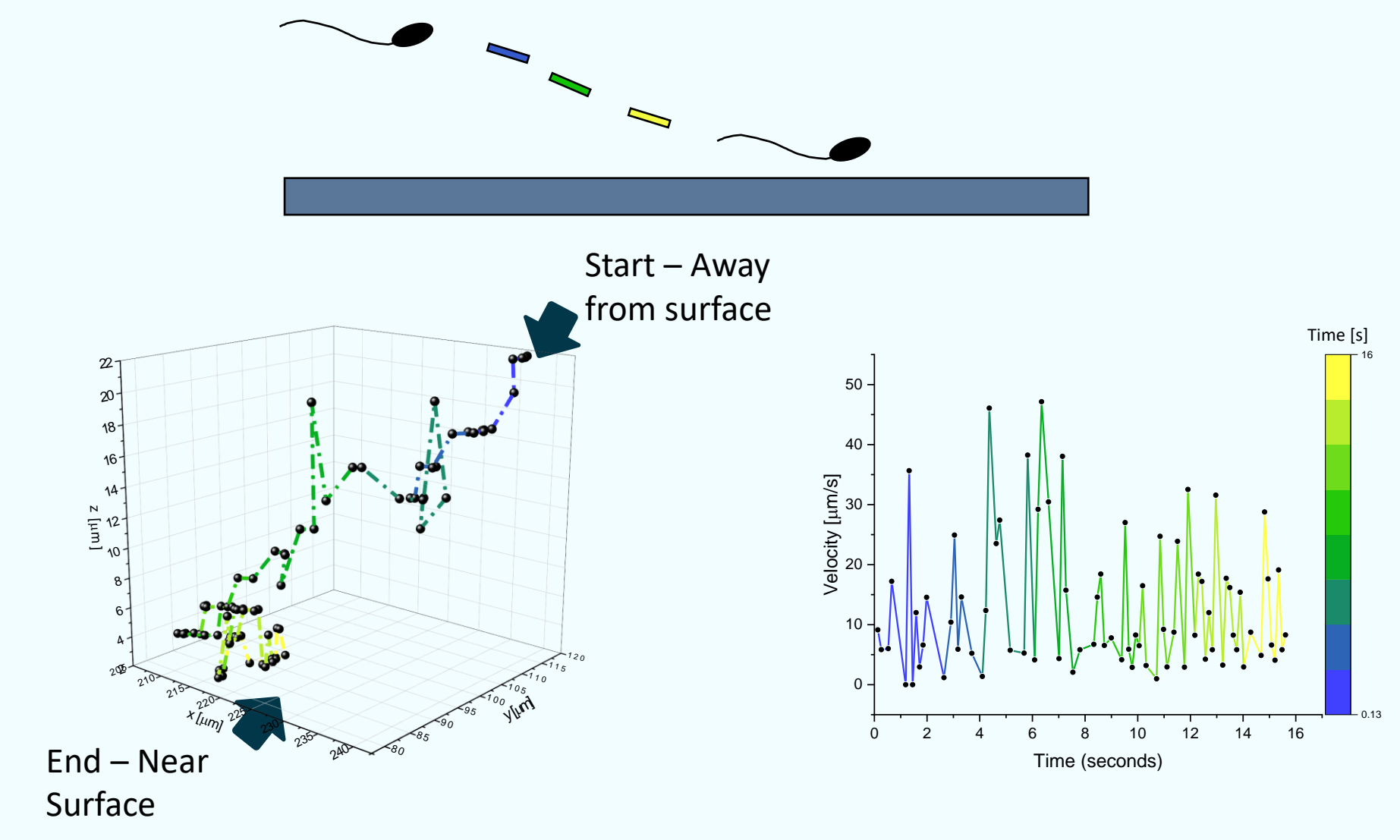

• Swimming speed transitions to lower speeds when swimming next to the surface.

4. Comparison of swimming speeds between normal gravity with rotation and simulated microgravity

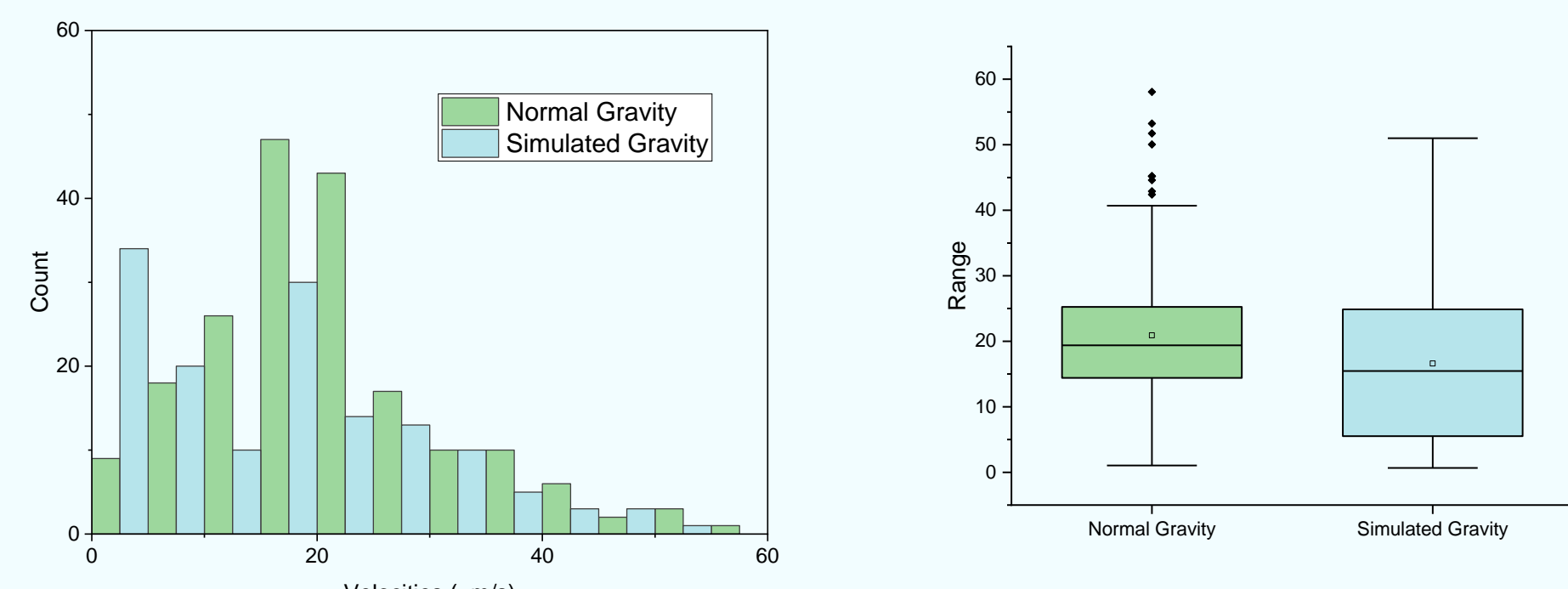

• Within error margins, these plots indicate little difference in swimming speeds between normal gravity and simulated microgravity

Acknowledgements

Dr. Jay Nadeau
Iulia Hanczareck
Gary Cadima (graphics)
JPL – OWLS
Oregon Space Grant Consortium (OSGC)